



ESSENTIAL EAR, NOSE AND THROAT INSTRUMENTS



**SHEFFIELD SURGICAL
INSTRUMENTS**

GENERAL INSTRUMENTS

Page



SCALPEL HANDLES _____	5
SCISSORS _____	5
NEEDLE HOLDERS _____	5
DISSECTING FORCEPS _____	6
ARTERY FORCEPS _____	6
SUTURE CARRIERS _____	7
TROCARS _____	7
SUCTION TUBES _____	7

AURAL (EAR) INSTRUMENTS



HOOKS _____	8
MYRINGOTOMES _____	8
PROBES _____	8
SPECULA _____	8
CROCODILE FORCEPS _____	9
PICKS _____	9

NASAL INSTRUMENTS



PUNCHES _____	10
DRESSING FORCEPS _____	10
ELEVATORS _____	10
CARTILAGE KNIVES _____	11
SNARES _____	11
SPECULA _____	12,13
TURBINATE FORCEPS _____	13

ORAL (MOUTH) INSTRUMENTS



MOUTH GAGS _____	14
BIPODS & PLATES _____	15
TONSIL DISSECTORS _____	15
TONSIL HOLDING FORCEPS _____	15
ADENOID CURETTES _____	16
TONSIL SNARES _____	16
GOUGES _____	17
MALLETS _____	17
SYRINGES _____	17
DIATHERMY NEEDLES & CABLES _____	18
BIPOLAR FORCEPS & CABLES _____	18



SHEFFIELD SURGICAL ESSENTIAL EAR, NOSE AND THROAT (E.N.T.) INSTRUMENTS



This catalogue has been produced to show some of the most popular items from our extensive range of E.N.T. Surgical Instruments.

By working closely with some of the worlds's leading E.N.T. consultant surgeons, we have selected a range of instruments to cover the most common procedures and applications for ear, nose and throat surgery.

We have compiled three popular sets as shown on page 3 and page 4. The contents of each set will give a representation of the main instruments used for these procedures. Please note, however, that this is only a guideline and practitioners may wish to choose whichever instruments suit their own clinical preference.

This is a small selection from our large portfolio of E.N.T. instruments. We also produce a wide selection of instruments for general and middle ear applications. Details on our complete range is available on request.

All our E.N.T. instruments are manufactured from the highest quality materials and to the most exacting tolerances to provide excellent clinical performance and years of trouble free service. All instruments are fully autoclaveable.

Our stainless steel finish instruments are provided with a matt or satin finish to reduce glare under strong illumination. We also provide a number of instruments (especially for middle ear procedures) with a black coated finish to eliminate reflected light completely. These "blackened" instruments are very popular for use under illuminated magnification.

Some instruments such as the Thudicum specula and the Abbey needle are available with a special insulated coating for diathermy or cautery procedures.

As with all SHEFFIELD reusable instruments our E.N.T. range is covered by our unrivalled "Lifetime Guarantee" against faulty manufacture and materials. If the product fails during its working life, due to a manufacturing or material defect, we will repair or replace it - free of charge, with the minimum of fuss or delay.

If you have not yet experienced the benefits of SHEFFIELD instruments, simply contact Sheffield Surgical for samples for evaluation, you will not be disappointed.

SHEFFIELD SURGICAL INSTRUMENTS

RECOMMENDED E.N.T. INSTRUMENT SETS

In conjunction with some of the U.K.'s leading E.N.T. consultants, we have prepared a small selection of recommended sets for use in common E.N.T. procedures. These are by no means an exhaustive list of either the ear, nose and throat instruments that we produce or the sets we can provide. They should however, provide a guideline for common procedures.

TONSILLECTOMY & ADENOIDECTOMY TRAY

	Page
21.860.10 Yankauer Suction Tube Complete	7
30.345.20 Doughty Tongue Plate Slotted 2.5" (6.25cm)	14
30.345.30 Doughty Tongue Plate Slotted 3" (7.5cm)	14
30.345.40 Doughty Tongue Plate Slotted 3.5" (9cm)	14
30.345.50 Doughty Tongue Plate Slotted 4" (10cm)	14
30.345.60 Doughty Tongue Plate Slotted 4.5" (11.5cm)	14
30.300.10 Davis Boyle Mouth Gag Frame Child	14
30.300.30 Davis Boyle Mouth Gag Frame Adult	14
30.430.05 Draffin Bipod 4 Rings 15" (Paediatric) Pair	15
30.430.10 Draffin Bipod 4 Rings 19" Pair	15
30.960.30 St. Clair Thomson Adenoid Curette Size 1.5 12mm x 8.5" (21.5cm)	16
30.960.40 St. Clair Thomson Adenoid Curette Size 2 14mm x 8.5" (21.5cm)	16
30.960.50 St. Clair Thomson Adenoid Curette Size 3 16mm x 8.5" (21.5cm)	16
30.680.20 Denis Browne Tonsil Hold Forceps Large 8" (20cm)	15
87.660.20 Bipolar Forceps Straight 7" (18cm) Broad 1.00mm	18
87.850.20 Bipolar Cable 3 Meter Euro Silicon	18
12.965.20 Mcindoe Scissors Curved 7" (18cm) Blunt/Blunt Tungsten Carbide	5
14.020.10 Birkett Artery Forceps Straight 7.5" (19cm)	6
14.380.50 Negus Artery Forceps Long Curved 7.5" (19cm)	6
30.580.10 Mollison Tonsil Retractor/Dissector 8.25" (21cm)	15
28.400.10 Hill Elevator & Dissector Double Ended 9" (23cm)	10
17.960.10 Negus Knot Tier And Pusher 8" (20cm)	7
13.425.30 Waugh Dissecting Forceps 8" (20cm) Serrated	6
13.425.60 Waugh Dissecting Forceps 8" (20cm) 1x2 Teeth Serrated	6

INSERTION OF GROMMET & EXAMINATION UNDER ANAESTHETIC (E.U.A.) TRAY

26.813.15 Gruber Speculum 3mm Diamond Oval Black	8
26.813.20 Gruber Speculum 4mm Diamond Oval Black	8
26.813.25 Gruber Speculum 5mm Diamond Oval Black	8
26.889.35 Tumarkin Speculum No:4 Round Black Slot	8
26.675.20 Jobson Horne Probe Double Ended 7" (18cm)	8
21.870.10 Zoellner Suction Tube	7
21.870.30 Zoellner Suction Tube End 18swg	7
27.240.40 Hartman Crocodile Forceps Micro Serrated Black	9
26.130.10 Formby Cerumen Scoop & Hook 7" (18cm)	8
27.320.30 Ormerod Crocodile Forceps Micro Cup Black	9
27.685.50 Schuknecht Pick Straight 178mm Black	9
26.560.10 Agnew (Poltzer) Myringotome 7" (18cm)	8



28.715.40	Thudicum Nasal Speculum Set of 5	12
28.685.10	Killian Nasal Speculum 2" (5cm) x 5.5" (13cm) Self Retaining	12
28.685.20	Killian Nasal Speculum 2.5" (6.5cm) x 5.5" (13cm) Self Retaining	12
28.685.30	Killian Nasal Speculum 3" (7.5cm) x 5.5" (13cm) Self Retaining	12
28.685.40	Killian Nasal Speculum 3.5" (9cm) x 5.5" (13cm) Self Retaining	12
28.750.20	Thudicum Nasal Speculum No:2 Diathermy Insulated	13
11.550.10	Scalpel Handle No:3 Stainless Steel	5
28.470.80	Freer Elevator Sharp & Dull Blade 7.5" (19cm) Double Ended (8)	10
28.400.10	Hill Elevator/Dissector Double Ended 9" (23cm)	10
28.450.10	Ballenger Nasal Knife Straight 6mm Blade 8" (20cm)	11
28.310.10	Tilley Henckel Punch 4.5mm Jaw x 4.5" (11.5cm) Shaft	10
28.310.20	Tilley Henckel Punch 5.5mm Jaw x 4.5" (11.5cm) Shaft	10
40.798.10	Small Mallet Stainless Steel 8.5oz Smooth Head	17
40.530.20	Jenkin Gouge Straight 6mm x 5.5" (14cm)	17
40.530.25	Jenkin Gouge Straight 8mm x 5.5" (14cm)	17
28.805.10	Luc Turbinate Forceps Oval Small 6mm x 7.5" (19cm)	13
28.805.20	Luc Turbinate Forceps Oval Medium 7.5mm x 7.5" (19cm)	13
28.805.30	Luc Turbinate Forceps Oval Large 9mm x 7.5" (19cm)	13
28.380.10	Tilley Nasal Dressing Forceps 7" (18cm)	10
21.870.10	Zoellner Suction Tube	7
21.870.35	Zoellner Suction Tube End 20swg	7
21.780.30	Self Clearing Suction Tube Medium	7
12.415.25	Heymann Turbinectomy Scissors Angled 7" (18cm)	5
12.985.20	Strabismus Scissors Curved 4.5" (11.5cm) Blunt/Blunt Tungsten Carbide	5
17.330.10	Crilewood Needle Holder 6" (15.5cm) Tungsten Carbide	5
14.325.20	Halstead Mosquito Artery Forceps Straight 5" (13cm) Standard	6
745.550.40	Dental Cartridge (A) Syringe Self Aspirating 2.2ml	17
86.800.40	Abbey Sub-Mucosal Needle Diathermy Insulated	18
86.251.20	Diathermy Cable 3 Meter Pin Fitting 4mm End	18
21.200.25	Tilley Litchwitz Trocar Complete 7.5fg	7
21.200.65	Tilley Litchwitz Cannula Only 7.5fg	7
28.620.10	Glegg Nasal Snare Self Forming 9" (23cm)	11
28.635.10	Negus Nasal Snare 9" (23cm)	11
13.225.10	Gillies Dissecting Forceps 6" (15.5cm) 1x2 Teeth Serrated	6
13.295.10	Mcindoe Dissecting Forceps 6" (15.5cm) Serrated	6

SCALPEL HANDLES



11.550.10 SCALPEL HANDLE NO.3 STAINLESS STEEL

11.550.30 SCALPEL HANDLE NO.4 STAINLESS STEEL (NOT ILLUSTRATED)

SCISSORS

Although we produce other Turbinectomy scissors, such as Beckman and Cottle, the Heyman is the most popular for Turbinectomy procedures.

Our full range of T-Carb Tungsten Carbide instruments are available in our full colour 24 page Tungsten Carbide Catalogue.



12.415.25 HEYMAN TURBINECTOMY SCISSORS ANGLED 7" (18CM)




12.965.20 MCINDOE SCISSORS CURVED 7" (18CM) BLUNT/BLUNT T.C.



12.985.20 STRABISMUS SCISSORS CURVED 4.5" (11.5CM) BLUNT/BLUNT T.C.

12.985.10 STRABISMUS SCISSORS STRAIGHT 4.5" (11.5CM) BLUNT/BLUNT T.C. (NOT ILLUSTRATED)

NEEDLEHOLDERS

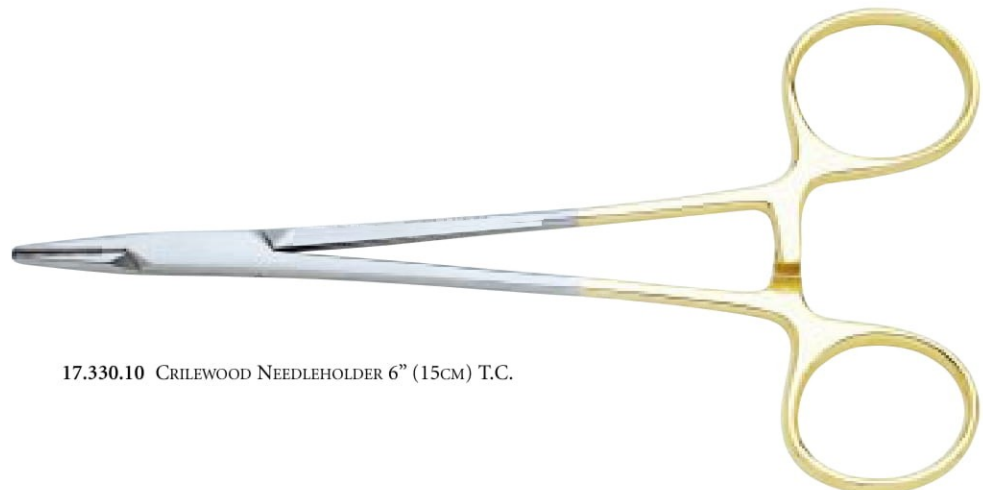

17.330.10 CRILEWOOD NEEDLEHOLDER 6" (15CM) T.C.



DISSECTING FORCEPS



13.225.10 GILLIES DISSECTING FORCEPS 6" (15.5CM) 1X2 TEETH SERRATED

13.226.10 GILLIES DISSECTING FORCEPS 6" (15.5CM) 1X2 TEETH SERRATED FINE (NOT ILLUSTRATED)



13.295.10 MCINDOE DISSECTING FORCEPS 6" (15.5CM) SERRATED

13.295.20 MCINDOE DISSECTING FORCEPS 6" (15.5CM) 1X2 TEETH (NOT ILLUSTRATED)



13.425.30 WAUGH DISSECTING FORCEPS 8" (20CM) SERRATED

(ALSO AVAILABLE BUT NOT ILLUSTRATED)

13.425.20 WAUGH DISSECTING FORCEPS 7" (18CM) SERRATED


13.425.50 WAUGH DISSECTING FORCEPS 7" (18CM) 1X2 TEETH SERRATED

13.425.60 WAUGH DISSECTING FORCEPS 8" (20CM) 1X2 SERRATED

ARTERY FORCEPS

14.020.10 BIRKETT ARTERY FORCEPS STRAIGHT 7.5" (19CM)

14.020.20 BIRKETT ARTERY FORCEPS CURVED 7.5" (19CM) (NOT ILLUSTRATED)

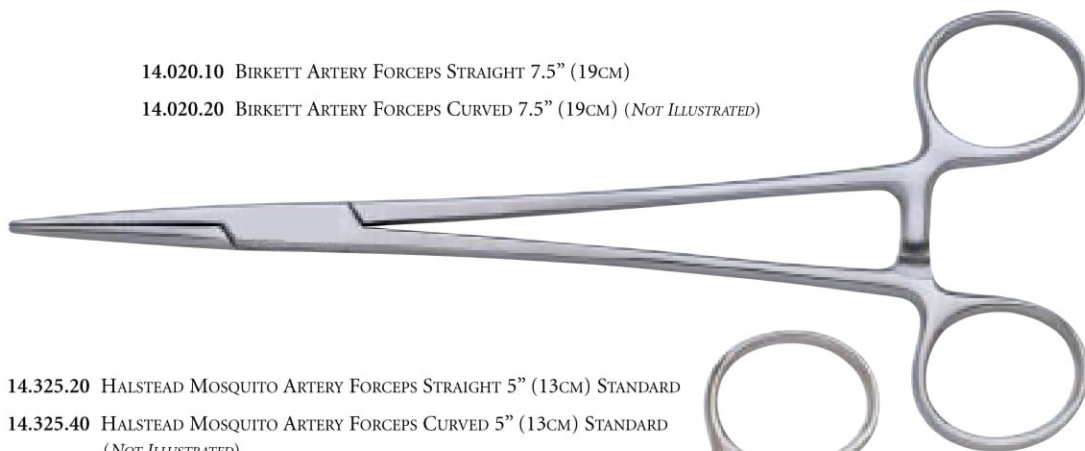

14.325.20 HALSTEAD MOSQUITO ARTERY FORCEPS STRAIGHT 5" (13CM) STANDARD

14.325.40 HALSTEAD MOSQUITO ARTERY FORCEPS CURVED 5" (13CM) STANDARD (NOT ILLUSTRATED)



14.380.50 NEGUS ARTERY FORCEPS LONG CURVED 7.5" (19CM)




14.380.40 NEGUS ARTERY FORCEPS SHORT CURVED 7.5" (19CM)

Negus - Two types of jaws are available for grasping the tonsil nodule. The long curve model is the most popular of the two versions.

SUTURE CARRIERS

Our Negus Knot Tier has a fine tip with smooth edges to avoid sutures snagging or cutting.

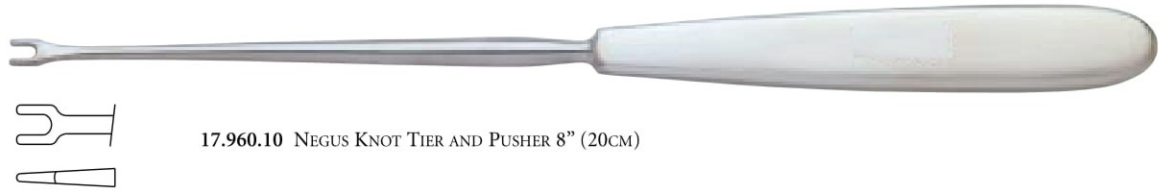

17.960.10 NEGUS KNOT TIER AND PUSHER 8" (20CM)

TROCARS

Traditionally used for antral wash outs. Other sizes are also available.



21.200.25 TILLEY LITCHWITZ TROCAR COMPLETE 7.5FG

21.200.20 TILLEY LITCHWITZ TROCAR COMPLETE 7FG (NOT ILLUSTRATED)

SUCTION TUBES

Self Clearing Suction Tube - This model has a built in stillette which allows the tube to be cleared of any blockages whilst in use.



21.780.20 - Small

21.780.30 - Medium

21.780.40 - Large

21.780.30 SELF CLEARING SUCTION TUBE MEDIUM

21.780.20 SELF CLEARING SUCTION TUBE SMALL

21.780.40 SELF CLEARING SUCTION TUBE LARGE

The Yankauer Suction Tube has a rose end which may be removed for cleaning. Ideal for suction in tonsil and adenoid procedures.



21.860.10 YANKAUER SUCTION TUBE COMPLETE

The Zoellner Suction Tube is ideal for small, localised suction due to its narrow size.

The various detachable suction tube ends allow

21.870.10

21.870.10 ZOELLNER SUCTION TUBE

21.870.30 ZOELLNER SUCTION TUBE END 18SWG

21.870.35 ZOELLNER SUCTION TUBE END 20SWG

21.870.40 ZOELLNER SUCTION TUBE END 22SWG

21.870.45 ZOELLNER SUCTION TUBE END 24SWG

21.870.50 ZOELLNER SUCTION TUBE END 26SWG



MYRINGOTOMES, HOOKS AND PROBES

The super sharp trocar point of the Agnew myringotome is ideal for easily perforating the ear drum.



26.560.10 AGNEW (POLITZER) MYRINGOTOME 7" (18CM)



26.130.10 FORMBY CERUMEN SCOOP + HOOK 7" (18CM)



26.675.20 JOBSON HORNE PROBE DOUBLE ENDED 7" (18CM)

SPECULA




26.813.10 GRUBER SPECULUM SET OF 4 OVAL BLACK

- 26.813.13 GRUBER SPECULUM 2MM OVAL BLACK
- 26.813.15 GRUBER SPECULUM 3MM OVAL BLACK
- 26.813.20 GRUBER SPECULUM 4MM OVAL BLACK
- 26.813.25 GRUBER SPECULUM 5MM OVAL BLACK
- 26.813.30 GRUBER SPECULUM 7MM OVAL BLACK



26.889.10 TUMARKIN SPECULUM SET OF 8, 0-7 ROUND BLACK SLOT

- 26.889.15 TUMARKIN SPECULUM No.0 ROUND BLACK SLOT
- 26.889.20 TUMARKIN SPECULUM No.1 ROUND BLACK SLOT
- 26.889.25 TUMARKIN SPECULUM No.2 ROUND BLACK SLOT
- 26.889.30 TUMARKIN SPECULUM No.3 ROUND BLACK SLOT
- 26.889.35 TUMARKIN SPECULUM No.4 ROUND BLACK SLOT
- 26.889.40 TUMARKIN SPECULUM No.5 ROUND BLACK SLOT
- 26.889.45 TUMARKIN SPECULUM No.6 ROUND BLACK SLOT
- 26.889.50 TUMARKIN SPECULUM No.7 ROUND BLACK SLOT

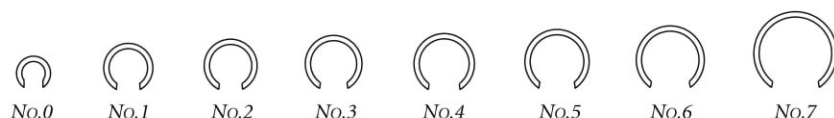

The Gruber and Tumarkin are our most popular ear specula. We have shown here the blackened variety, which are specially coated to prevent reflection or glare experienced when using the non coated steel versions.

FORCEPS AND PICKS



We manufacture a large range of crocodile forceps for middle ear procedures. Various models are available for grasping, cutting, biopsies and grommet insertions. A small selection of the blackened versions are shown here, which reduce glare or reflection under bright illumination. We also supply these in uncoated steel versions.



29

27

26

25

24

27.240.40 HARTMAN CROCODILE MICRO SERRATED BLACK

23

22

27.245.40 HARTMAN CROCODILE SUPER FINE SERRATED BLACK

21

20

27.246.40 HARTMAN CROCODILE EXTRA SUPER FINE SERRATED BLACK

19

18

17

16

27.320.30 ORMEROD CROCODILE MICRO CUP BLACK

15

14

27.320.20 ORMEROD CROCODILE MICRO SERRATED UP BLACK

13

12

27.320.40 ORMEROD CROCODILE MICRO ELONGATED CUP BLACK

11

10

27.320.50 ORMEROD CROCODILE MICRO SCISSORS BLACK

9

8

7

6

5

27.685.50 SCHUKNECHT PICK STRAIGHT 178MM BLACK

4

3


2

9



PUNCHES AND DRESSING FORCEPS

28.380.10 TILLEY NASAL DRESSING FORCEPS 7" (18CM)

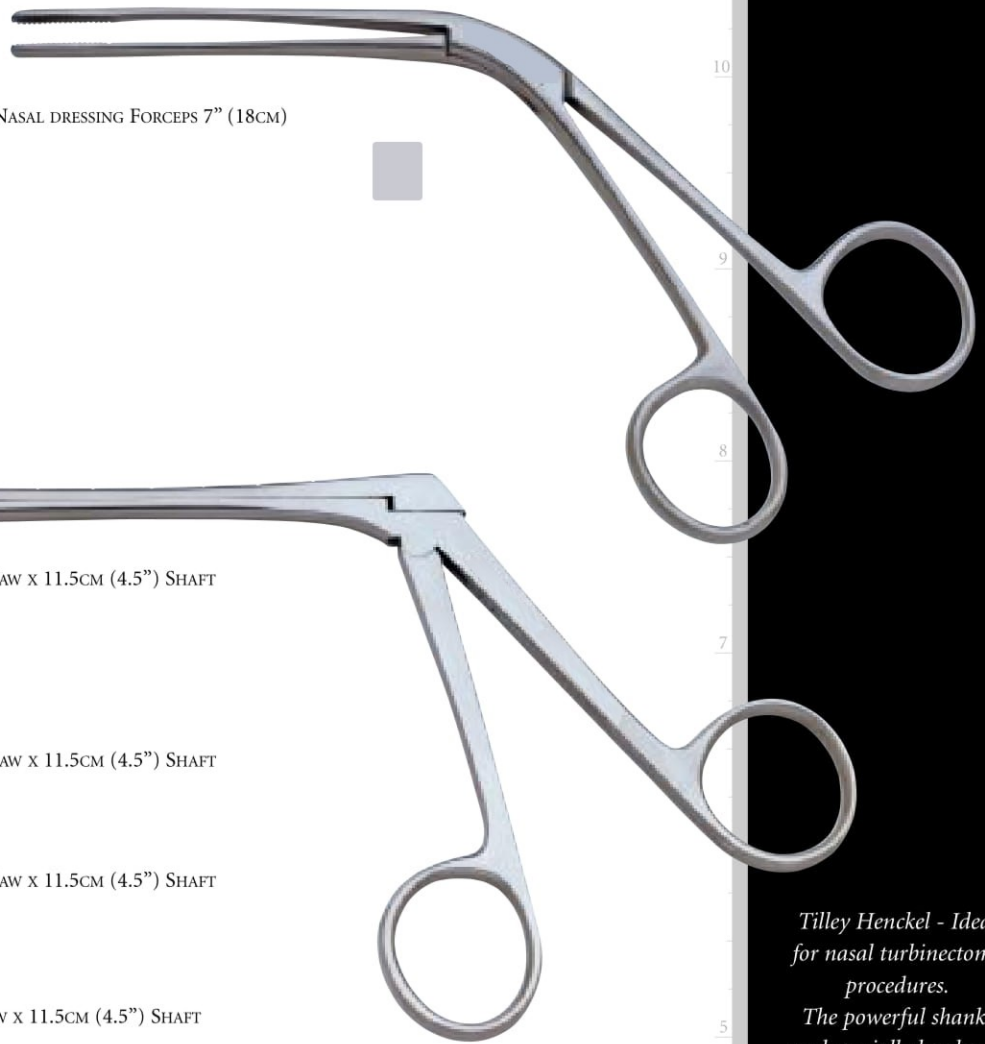

28.310.10 TILLEY HENCKEL PUNCH 4.5MM JAW X 11.5CM (4.5") SHAFT



28.310.20 TILLEY HENCKEL PUNCH 5.5MM JAW X 11.5CM (4.5") SHAFT



28.310.30 TILLEY HENCKEL PUNCH 6.5MM JAW X 11.5CM (4.5") SHAFT



28.310.40 TILLEY HENCKEL PUNCH 8MM JAW X 11.5CM (4.5") SHAFT

Tilley Henckel - Ideal for nasal turbinectomy procedures. The powerful shanks and specially hardened jaws provide a strong, clean bite.

ELEVATORS



28.400.10 HILL ELEVATOR/ DISSECTOR DOUBLE ENDED 9" (23CM)




28.470.80 FREER ELEVATOR SHARP + BLUNT BLADE 7.5" (19CM) DOUBLE ENDED

We manufacture a full range of Freer elevators, this double ended model is the most popular model with a sharp and a blunt blade.




CARTILAGE KNIVES


Ballenger Knives have super sharp blades to allow easy and controlled cutting of the nasal cartilage. The angled version (28.450.30) has a crank bend in the shank for use where access may be restricted.

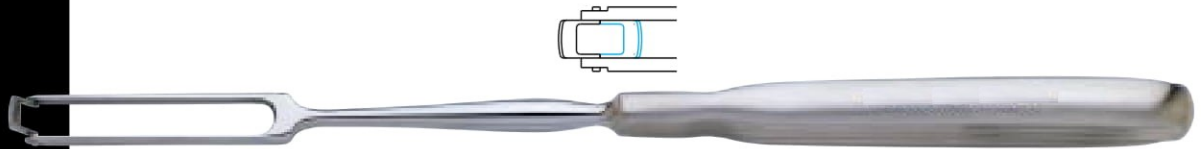


28.450.10 BALLENGER NASAL KNIFE STRAIGHT 8" (20CM) 6MM BLADE

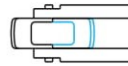

28.450.05 BALLENGER NASAL KNIFE STRAIGHT 8" (20CM) 4MM BLADE

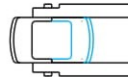

28.450.20 BALLENGER NASAL KNIFE STRAIGHT 8" (20CM) 8MM BLADE

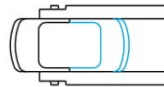

28.450.30 BALLENGER NASAL KNIFE ANGLED 8" (20CM) 6MM BLADE

NASAL SNARES



The Glegg Nasal Snare has a mechanism in the tip which pushes the wire apart as it is


extended from the tip, thus automatically forming a loop. The Negus requires the loop to be formed manually by hand before use. Both snares can be disassembled for replacement of the wire



28.620.10 GLEGG NASAL SNARE SELF FORMING 9" (23CM)



28.635.10 NEGUS NASAL SNARE 9" (23CM)

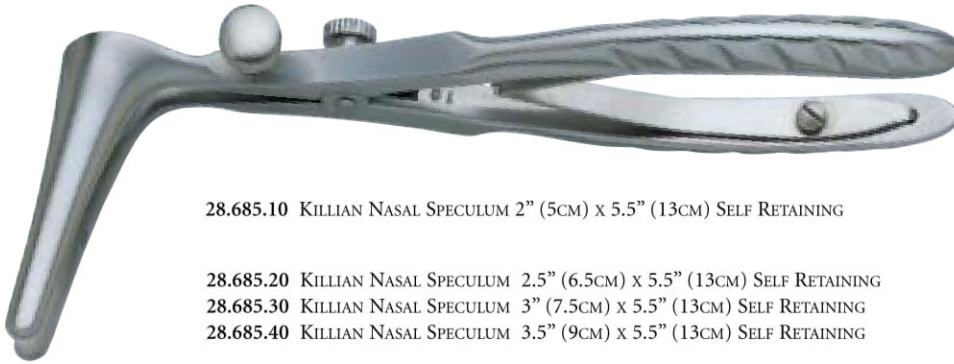


28.685.10 KILLIAN NASAL SPECULUM 2" (5CM) X 5.5" (13CM) SELF RETAINING

28.685.20 KILLIAN NASAL SPECULUM 2.5" (6.5CM) X 5.5" (13CM) SELF RETAINING

28.685.30 KILLIAN NASAL SPECULUM 3" (7.5CM) X 5.5" (13CM) SELF RETAINING

28.685.40 KILLIAN NASAL SPECULUM 3.5" (9CM) X 5.5" (13CM) SELF RETAINING

The Killian is a very popular deep nasal speculum due to the various lengths available and its self retaining mechanism. Its simple thumb screw, allows it to stay open, unaided whilst in use.

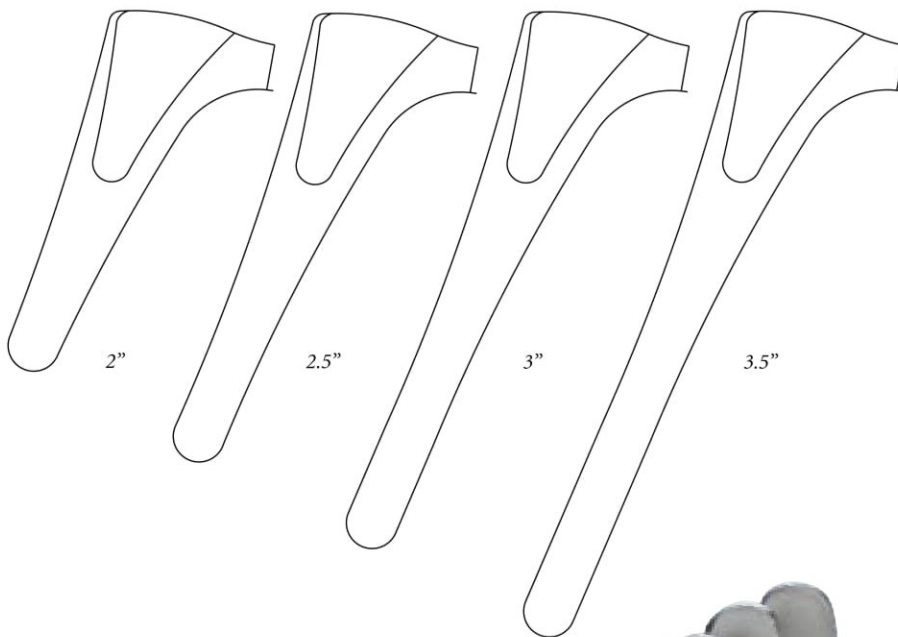

2"

2.5"

3"

3.5"

The Thudicum specula are ideal for general nasal procedures. They are available in a wide selection of sizes in both steel and an insulated diathermy coating.

28.715.40 THUDICUM NASAL SPECULUM SET OF 5




28.715.10 THUDICUM NASAL SPECULUM No.0

28.715.15 THUDICUM NASAL SPECULUM No.1

28.715.20 THUDICUM NASAL SPECULUM No.2

28.715.25 THUDICUM NASAL SPECULUM No.3

28.715.30 THUDICUM NASAL SPECULUM No.4

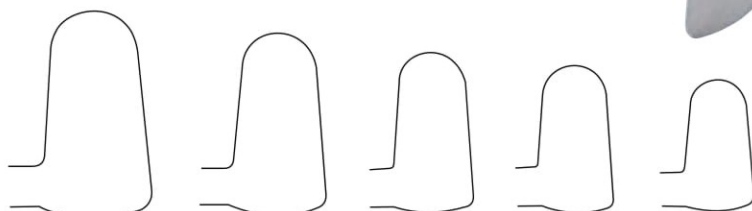

No.4

No.3

No.2

No.1


No.0

THUDICUM SIZES

NASAL SPECULAE

Our insulated diathermy Thudicums may be used with electrosurgery procedures. The insulated coating is fully autoclaveable.

The Mark Hovell variety are similar in design to the Thudicums but with a slightly different head shape and size. These are also available in insulated diathermy versions.

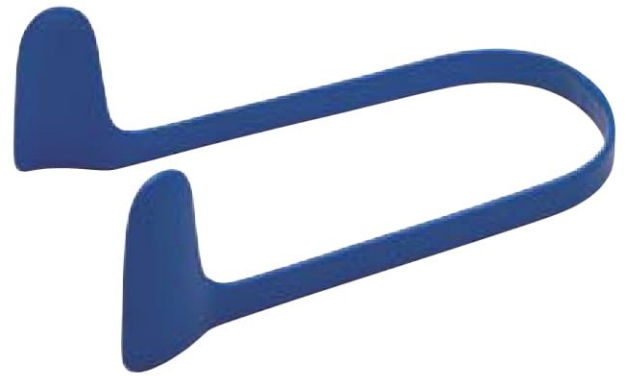


28.750.20 THUDICUM NASAL SPECULUM No.2 INSULATED

- 28.750.10 THUDICUM NASAL SPECULUM No.0 INSULATED
- 28.750.15 THUDICUM NASAL SPECULUM No.1 INSULATED
- 28.750.25 THUDICUM NASAL SPECULUM No.3 INSULATED
- 28.750.30 THUDICUM NASAL SPECULUM No.4 INSULATED
- 28.750.50 THUDICUM NASAL SPECULUM SET OF 5 INSULATED

- 28.690.10 MARK HOVELL NASAL SPECULUM No:1
- 28.690.15 MARK HOVELL NASAL SPECULUM No:2
- 28.690.20 MARK HOVELL NASAL SPECULUM No:3
- 28.690.25 MARK HOVELL NASAL SPECULUM No:4
- 28.690.30 MARK HOVELL NASAL SPECULUM No:5
- 28.690.35 MARK HOVELL NASAL SPECULUM No:6
- 28.690.40 MARK HOVELL NASAL SPECULUM No:7
- 28.690.50 MARK HOVELL NASAL SPECULUM SET OF 7

- 28.740.10 MARK HOVELL NASAL SPECULUM No.1 INSULATED
- 28.740.15 MARK HOVELL NASAL SPECULUM No.2 INSULATED
- 28.740.20 MARK HOVELL NASAL SPECULUM No.3 INSULATED
- 28.740.25 MARK HOVELL NASAL SPECULUM No.4 INSULATED
- 28.740.30 MARK HOVELL NASAL SPECULUM No.5 INSULATED
- 28.740.35 MARK HOVELL NASAL SPECULUM No.6 INSULATED
- 28.740.40 MARK HOVELL NASAL SPECULUM No.7 INSULATED
- 28.740.60 MARK HOVELL NASAL SPECULUM SET OF 7 INSULATED

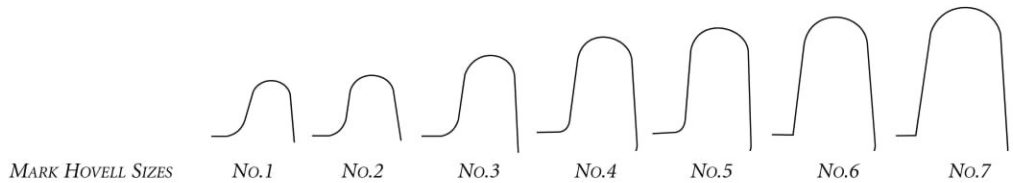


TURBINATE FORCEPS

Luc Turbinate - A very popular instrument for nasal turbinectomy procedures, the long jaws allow deep access and the super sharp jaws provide a clean and easy cut.



28.805.10 LUC TURBINATE FORCEPS OVAL SMALL 6MM X 7.5" (19CM)

- 28.805.20 LUC TURBINATE FORCEPS OVAL MEDIUM 7.5MM X 7.5" (19CM)
- 28.805.30 LUC TURBINATE FORCEPS OVAL LARGE 9MM X 7.5" (19CM)

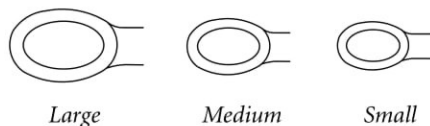

The Boyle Davis is the most popular mouth gag for Tonsillectomy and Adenoidectomy procedures. We supply a selection of tongue plates but the Doughty version shown here is the most commonly used. The slot sizes on our Doughty plates are made to the originally Doughty specifications.



30.300.10 DAVIS BOYLE MOUTH GAG FRAME CHILD

30.300.30 DAVIS BOYLE MOUTH GAG FRAME ADULT (NOT ILLUSTRATED)



WE ALSO SUPPLY SPARE TEETH GUARDS/PROTECTORS FOR MOUTH GAGS, PLEASE CONTACT US FOR DETAILS OR REFER TO OUR PRICE LIST


30.345.20 DOUGHTY TONGUE PLATE 2.5" (6.25CM) WITH 2MM SLOT

30.345.30 DOUGHTY TONGUE PLATE 3" (7.5CM) WITH 3MM SLOT

30.345.40 DOUGHTY TONGUE PLATE 3.5" (9CM) WITH 5.7MM SLOT

30.345.50 DOUGHTY TONGUE PLATE 4" (10CM) WITH 7MM SLOT

30.345.60 DOUGHTY TONGUE PLATE 4.5" (11.5CM) WITH 8MM SLOT

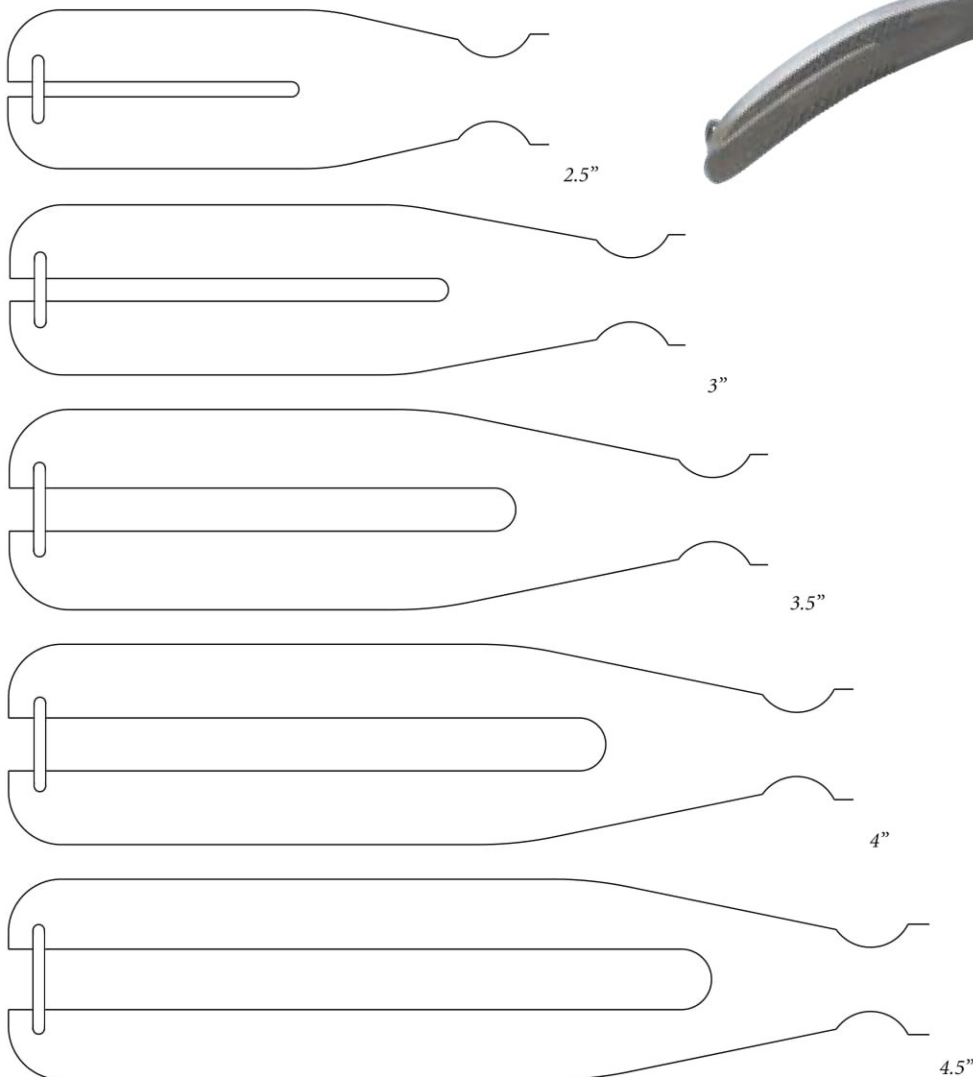

BIPODS AND PLATES

Our Draffin Bipods are made from specially hardened stainless steel to ensure they do not bend or deform whilst in use.



30.430.05 DRAFFIN BIPOD 4 RINGS 15" (38.5CM) PAIR

30.430.10 DRAFFIN BIPOD 4 RINGS 19" (48CM) PAIR (NOT ILLUSTRATED)



The Magauran plates are used with both child (15") and adult (19") bipods.

30.440.10 MAGAURAN PLATE 16" (40.5CM)

TONSIL DISSECTORS AND FORCEPS



30.580.10 MOLLISON TONSIL RETRACTOR/DISSECTOR 8.25" (21CM)


Our Denis Browne forceps have a fenestrated smooth cup to cradle the tonsil without perforating it.



30.680.20 DENIS BROWNE TONSIL HOLDING FORCEPS LARGE 8" (20CM)




ADENOID CURETTES

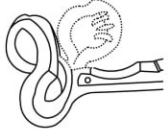

30.960.30 ST. CLAIR THOMSON ADENOID CURETTE 1.5 12MM x 8.5" (21.5CM)

- 30.960.10 ST. CLAIR THOMSON ADENOID CURETTE 0 8.5" (21.5CM) x 8MM
- 30.960.20 ST. CLAIR THOMSON ADENOID CURETTE 1 8.5" (21.5CM) x 10MM
- 30.960.40 ST. CLAIR THOMSON ADENOID CURETTE 2 8.5" (21.5CM) x 14MM
- 30.960.50 ST. CLAIR THOMSON ADENOID CURETTE 3 8.5" (21.5CM) x 16MM
- 30.960.60 ST. CLAIR THOMSON ADENOID CURETTE 4 8.5" (21.5CM) x 18MM

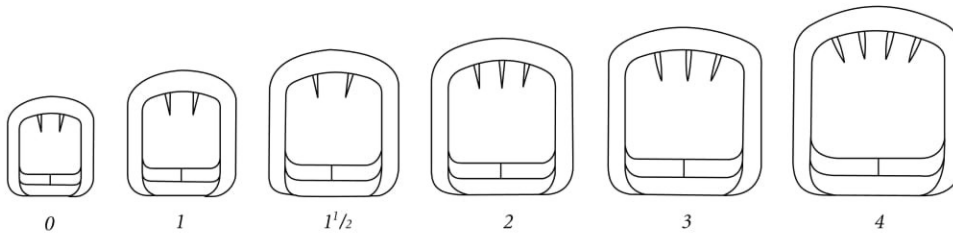

The Cage on the jaw of the St Clair Thomson Currettes features spikes to securely hold the adenoids whilst the cutting blade slices it. The cage assembly may be disassembled for cleaning and blade sharpening.

TONSIL SNARE



30.895.10 EVE TONSIL SNARE 11" (28CM)

The Eve Tonsil snare is ideal for tonsillectomies. The snare can be easily disassembled for replacing the snare wire.



GOUGES AND MALLETS

The Jenkins gouge is suited for nasal bone fractures and removal of the nasal hump.

Other sizes are also available, as are other nasal gouges such as Tilley, Killian and Guy.



40.530.20 JENKIN GOUGE STRAIGHT 6MM X 5.5" (14CM)

40.530.25 JENKIN GOUGE STRAIGHT 8MM X 5.5" (14CM)

This mallet is also commonly referred to as a Toffee Hammer or Mallet. Due to its light weight and small size it is easy to control and is popularly used with chisels, osteotomes and gouges for a variety of minor orthopaedic procedures.




40.798.10 SMALL MALLET STAINLESS STEEL 8.5OZ SMOOTH HEAD

SYRINGES

This syringe is designed to be used with rubber stopped cartridges.


The syringe will self aspirate when used with the correct cartridges.

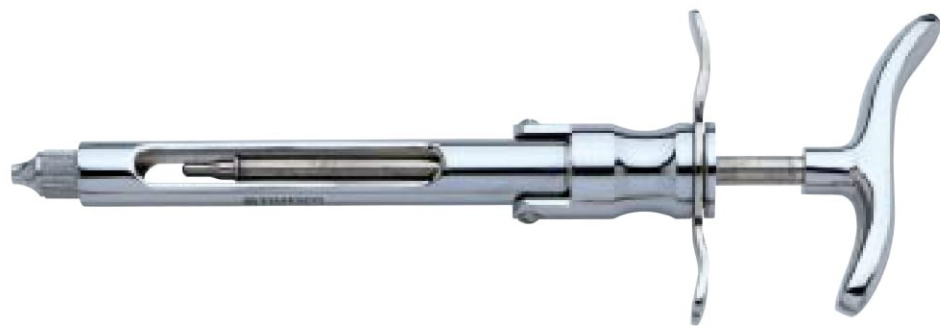

745.550.40 DENTAL CARTRIDGE (A) SYRINGE SELF ASPIRATING 2.2ML



DIATHERMY NEEDLES & CABLES

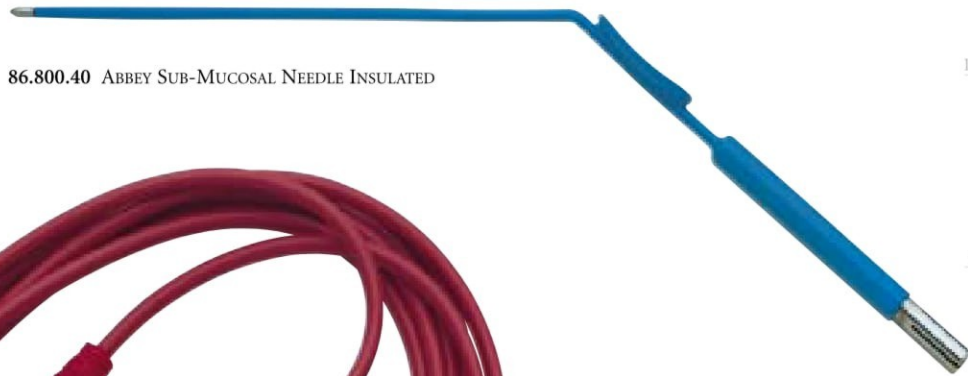

86.800.40 ABBEY SUB-MUCOSAL NEEDLE INSULATED



86.251.20 DIATHERMY CABLE 3 METRES PIN (4MM) END

86.251.40 DIATHERMY CABLE 3 METRES PIN (8MM) END (NOT ILLUSTRATED)

The Abbey Needle has an insulated diathermy coating to allow it to be used for mucosal cauterisation.

Our diathermy cable is designed to be used with diathermy or monopolar forceps and instruments. Available in 2 pin sizes, 4mm or 8mm.

BIPOLAR FORCEPS & CABLES



87.660.20 BIPOLAR FORCEPS STRAIGHT 7" (18CM) BROAD 1.00MM TIP

87.850.20 BIPOLAR CABLE 3 METRES EURO SILICON




A full range of Bipolar Forceps is available in our Bipolar Forceps brochure.

Our Bipolar cable is designed for use with bipolar forceps and instruments. Available with a European fitting (as shown) or a U.S.A. two pin fitting (not shown).





SURGICAL

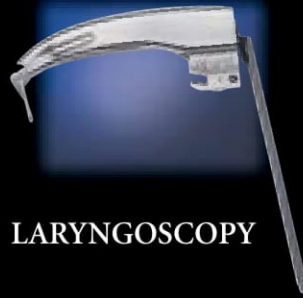

LARYNGOSCOPY



PRE-HOSPITAL



FOOTCARE



DENTAL



**SHEFFIELD
SURGICAL
INSTRUMENTS**

Quality

Instruments

Since 1992

Sheffield Surgical Building, Wazirabad Road, Ugoki, Sialkot-51310 Pakistan.

Tell # +92-333-1400682

Email: sheffieldsurgicalinstruments@gmail.com

Email: info@sheffieldsurgicalinst.com

Web: www.sheffieldsurgicalinst.com

